



# CR 85-X

## Digitizer

*CR 85-X is a multi-user digitizer featuring a unique drop-and-go buffer that eliminates waiting times and maximizes productivity*

*CR 85-X is a multi-application digitizer, benefiting from three different image resolution modes*

### MAXIMIZING PRODUCTIVITY FOR THE COMPLETE RANGE OF CLINICAL APPLICATIONS

- No waiting times for improved patient care
- Input/output buffer for maximized productivity
- For a broad range of applications
- Three different image resolution modes (pixel pitch [ $\mu\text{m}$ ] : 50 - 100 - 150)

#### Highest productivity

The cassette buffer eliminates waiting times and allows for a continuous workflow within the department. Zero-button operation with automated cassette handling makes CR 85-X a highly productive and user-friendly system with a throughput of up to 115 plates an hour, depending on size and application. Using CR 85-X as a central digitizer in the radiology department, multiple examination rooms can be supported.

#### No waiting

The CR 85-X digitizer requires no manual interaction and all the user has to do is to deposit the cassettes in the input buffer (up to 10 cassettes). The digitizer automatically takes cassettes from the input buffer and reads the demographic data from the memory on the cassette. It then scans the imaging plate, digitizes the image and returns the cassette to the output buffer for new exposures.



*Integrated CR user station for time-saving identification and optimized workflow*

## Full data

CR 85-X reads imaging plates at a standard resolution of 6 pixels/mm. The high resolution mode of 10 pixels/mm is available for all image plate sizes.

The maximum resolution mode of 20 pixels/mm is available for dedicated 18 x 24 cm and 24 x 30 cm extremities cassettes and plates.

## Broad range of applications

In combination with the application-specific plates and cassettes, CR 85-X supports a broad range of applications :

- General radiography
- Orthopaedics – Extremities
- Dental
- Paediatrics

## Compact footprint & optimal accessibility

CR 85-X occupies a very small floorspace and at the same time provides unhindered access to several users, both at the input and the output buffer, resulting in a smooth flow of operations. This concept makes CR 85-X the state-of-the-art solution for centralized CR environments.



## CR user station

Its modular and ergonomic design includes:

- Cassette identification functions
- Space for:
  - Workstation for image handling, processing and dispatching
  - Monitor, network switches and UPS
  - Cassette storage.


## An economical way to go digital

CR is compatible with all existing X-ray systems allowing X-ray departments to go digital without significant additional investments and workflow adaptations.



## CASSETTE SIZES: CR 85-X

Accepted Cassette Sizes	Spatial Resolution	Pixel Matrix Size
<b>Standard resolution</b>		
35 x 43 cm (14 x 17 in)	6 pixels / mm	2320 x 2828
35 x 35 cm (14 x 14 in)	6 pixels / mm	2320 x 2320
<b>High resolution</b>		
35 x 43 cm (14 x 17 in)	10 pixels / mm (option)	3480 x 4240
35 x 35 cm (14 x 14 in)	10 pixels / mm (option)	3480 x 3480
35 x 43 cm (automatic collimation to 21 x 43 cm)	10 pixels / mm	2020 x 4240
24 x 30 cm	10 pixels / mm	2320 x 2920
18 x 24 cm	10 pixels / mm	1720 x 2320
15 x 30 cm	10 pixels / mm	1420 x 2920
8 x 10 in	10 pixels / mm	1950 x 2460
10 x 12 in	10 pixels / mm	2460 x 2970

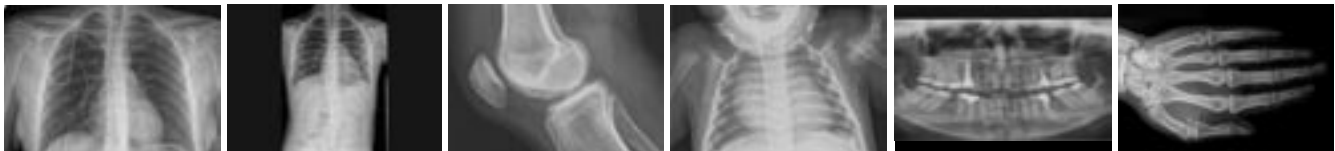

## SAFETY

Region	Regulation	X-ray	Laser
<b>Europe</b>	EN 60601-1: 1990 + A1: 1993 + A2: 1995 EN 60601-1-2: 2001	Regulation: 1987	EN 60825 - 1:2001
<b>USA</b>	UL 2601 21CFR part 820: good manufacturing practice for medical devices	DHHS/FDA 21 CFR part 1002, subchapter B	DHHS/FDA 21 CFR parts 1040, 10 and 1040, 11
<b>Canada</b>	CSA22.2 No.601.1 No.601.1.2		

# technical

## SPECIFICATIONS

### GENERAL

#### Cassette buffer capacity and performance

- 10 cassettes of mixed sizes, both in input and output buffer
- throughput: up to 112 plates/h (depending on size and application)

#### LCD display

- Machine status and error conditions

#### Greyscale resolution

- Data acquisition: 12 bits/pixel
- Output to processor: 12 bits/pixel

#### Dimensions and weight

- W x d x H: 84 x 115 x 142 cm (33 x 45 x 56 in)
- At foot: 84 cm (33 in)
- At buffer: 142 cm (56 in)
- Weight: Approx. 397 kg (875,2 lbs)

#### Power

- 50/60 Hz single phase
- 240V +10%, max. fuse 16A

- 230V ±10%, max. fuse 16A
- 208V ±10%, max. fuse 15A (e.g. USA)
- 200V ±10%, max. fuse 15A (e.g. Japan)

#### Environmental conditions

- Temperature: 15 - 30 °C (59 - 86°F)
- Humidity: 15 - 80% RH
- Magnetic fields: max. 12.60 μT
- Rate of change of temperature: 0.5°C/minute (0.9°F)

#### Environmental effects

- Noise level: max. 65 dB (A)
- Heat dissipation: standby 350 W, max. 2000 W

### SAFETY

#### Approvals

- TÜV, UL, cUL, CE

#### Transport details

- Temperature: -25 to +55°C (-4 to 131°F), -25°C for max. 72 hours, +55°C for max. 96 hours
- Humidity: 5 - 95% RH



Agfa-Gevaert has been certified by Lloyd's Register Quality Assurance Limited to the following quality management system standard: ISO 9001:2000. The quality management system is applicable to: For Healthcare applications - Marketing, design, development and production of imaging and communication solutions (film, paper and plates, chemicals, components, equipment and software).



Agfa-Gevaert has been awarded the Approval of Conformity certificate by Lloyd's Register Quality Assurance. It certifies that the Quality Management System for our X-Ray films conforms to the requirements of Annex V of the EEC Directive 93/42. Agfa HealthCare has been certified by Lloyd's Register Quality Assurance Limited to the following quality management system standards: ISO 9001:2000 and EN ISO 13485-2003. The quality management system is applicable to: marketing, design, development and production of imaging and communication solutions (film, paper and plates, chemicals, components, equipment and software) for Healthcare applications. Agfa-Gevaert had been awarded the ISO 9001 certificate by TÜV Zertifizierungsgemeinschaft e.V. This is applicable to Agfa's Quality Management System for design, production and servicing of Agfa Medical Equipment. Products distributed in North America are manufactured by/for Agfa Corporation, 10 South Academy St., Greenville, South Carolina 29601. Agfa, the Agfa rhombus and MUSICA are trademarks of Agfa-Gevaert N.V., Belgium or its affiliates. All other trademarks are held by their respective owners and are used in an editorial fashion with no intention of infringement. The data in this publication are for illustration purposes only and do not necessarily represent standards or specifications which must be met by Agfa. All information contained herein is intended for guidance purposes only, and characteristics of the products described in this publication can be changed at any time without notice. Products may not be available for your local area. Please contact your local sales representative for availability information. Agfa diligently strives to provide as accurate information as possible, but shall not be responsible for any typographical error.

### Agfa Corporation

HealthCare Headquarters  
10 South Academy Street,  
Greenville, SC 29601

### Canada

Agfa, Inc.  
77 Belfield Road, Etobicoke, Ont.,  
M9W 1G6

### Mexico

Agfa de México, S.A. de C.V.  
Benjamin Franklin #98, Colonia  
Escandon, C.P. 11800, México, D.F.