

Accommodates
BARIATRIC
Patients



BIODEX MEDICAL

NEW

Ultra Pro™ Ultrasound Table

Sonographer-friendly design.

- Ergonomically designed to help prevent musculoskeletal injuries
- 500-lb patient capacity
- Open frame design for greater accessibility
- Motor-controlled adjustments – including Fowler back, height, Trendelenburg motions and auto level
- Vascular scanning arm board (optional) can be mounted to either side of the table
- Side rails (optional) fold out of the way to allow unobstructed patient access
- Equipped with hand controller to actuate table positions
- Optional table accessories include foot controller, retractable stirrups, vascular scanning arm board, fold away side rails, IV pole and paper dispenser



The Ultra Pro Ultrasound Table is a dedicated table for all general ultrasound procedures including OB/GYN. Choose from an assortment of table options to customize the table to suit your department needs.

With a vast height range down to 23", the Ultra Pro easily accommodates wheelchair transfers,

patients with ambulatory difficulties and the bariatric population.

Ergonomically designed to suit both patient and sonographer, the table's height range allows seated or standing work positions while the fold away side rails allow comfortable shoulder, arm and hand positions when scanning. There are no barriers between sonographer and patient, allowing correct scanning posture to prevent musculoskeletal injury. The open frame design allows sonographers to get closer to the patient, especially useful for carotid imaging or long leg arterial scanning. The motorized Fowler back is infinitely adjustable up to 80° via hand or foot controller.

For the patient, comfort is assured through lengthy procedures with an amply thick mattress and a multitude of comfort-driven design details. For example, the retractable stirrups can extend extra long to accommodate taller or larger patients; the leg section drops down to 40° for patients with circulatory issues and 80° for stirrup access. Biodex tables are so accommodating, they actually reduce the total amount of time required to achieve a quality image.

The Ultra Pro Ultrasound Table is equipped with a motor-controlled hand controller to actuate table positions, 2"- and 3"-thick table padding covered in easy-to-clean Naugahyde®, two

patient restraining straps with Velcro® hook and loop fasteners, five-inch locking swivel casters and a central locking system, accessible from either side of the table, to secure table position.

Additional accessories are available, including foot controller, retractable stirrups, vascular scanning arm board, folding side rails and IV pole.

SPECIFICATIONS:

Dimensions: 70" l x 30" w (185.4 x 76.2 cm); 35" w (88.9 cm) with side rails;

Drop-Down Leg Section: 12.6" l (32 cm)

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Three primary sections: torso section with Fowler back positioning, center section remains fixed and leg section drops down to 40° and 80° for stirrup access

Mattress: Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) Individual locking swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat


Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 335 lb (152 kg)

Power: 115 VAC or 230 VAC

Warranty: One year

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2 No. 601.1-M90 and EN 60601-1 standards.

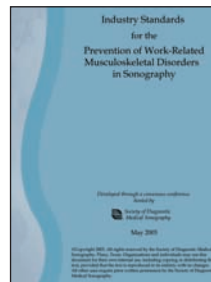

More than 80% of sonographers are scanning in pain and 20% of these professionals eventually experience a career-ending injury. On average, within five years of entering the profession, sonographers experience pain while scanning.

For more information on SDMS Industry Standards for the Prevention of Work-Related Musculoskeletal Disorders in Sonography, visit our website at www.biodex.com

058-720 Table, Ultrasound, Ultra Pro, 115 VAC
Includes Hand Controller.

*All options available separately**

058-725 Table, Ultrasound, Ultra Pro, 230 VAC
Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-632 Side Rails, Folding

058-651 Stirrups, Retractable

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-610 Dispenser, Paper
with Cutter Strap

Replacement:

056-612 Table Paper, Hygienic, 12 rolls
18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*