


Accommodates
BARIATRIC
Patients

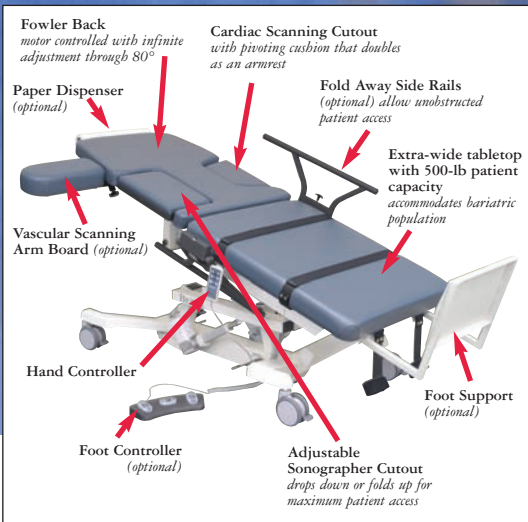

MEDICAL

NEW

Echo Pro™ Echocardiography Table

Efficient positioning reduces sonographer injuries and shortens exam time.

- Ergonomically designed to help prevent musculoskeletal injuries
- Extra wide top with 500-lb patient capacity – perfect for bariatric patients
- Motor-controlled adjustments – including Fowler back, height, Trendelenburg motions and auto level
- Cardiac scanning cutout features a patented swing-out cushion that doubles as an armrest
- Adjustable sonographer cutout allows comfortable scanning positioning
- Vascular scanning arm board (optional) can be mounted to either side of the table
- Side rails (optional) fold out of the way to allow unobstructed patient access
- Equipped with hand controller to actuate table positions
- Optional table accessories include foot controller, vascular scanning arm board, folding side rails, IV pole, paper dispenser and foot support



The Echo Pro Echocardiography Table is a dedicated solution for all echocardiography procedures. Extra wide with 500-lb patient capacity and a vast height range down to 23", this table will easily accommodate wheelchair transfers, patients with ambulatory difficulties and the bariatric population.

Ergonomically designed to suit both patient and sonographer, the

table's height range allows seated or standing work positions while the fold away side rails and adjustable sonographer cutout cushion allow comfortable shoulder, arm and hand positions when scanning. There are no barriers between sonographer and patient, allowing correct scanning posture to prevent musculoskeletal injury. The cardiac scanning cutout allows open access to the patient's left thorax area for an unobstructed apical approach. The motorized Fowler back is infinitely adjustable up to 80° via hand or foot controller. The table is so accommodating, stress-echo patients are off the treadmill, onto the table and into position for imaging quickly and comfortably.

For the patient, comfort is assured through lengthy procedures with cushioned arm supports, an amply thick mattress and a multitude of comfort-driven design details. For example, the unique swing-out cushion of the cardiac scanning cutout doubles as an armrest. The extra-wide tabletop comfortably accommodates the bariatric patient. In fact, the table is so accommodating, it reduces the total amount of time required to achieve a quality image.

The Echo Pro Echocardiography Table is equipped with a motor-

controlled hand controller to actuate table positions, 2"- and 3"-thick table padding covered in easy-to-clean Naugahyde®, two patient restraining straps with Velcro® hook and loop fasteners, five-inch locking swivel casters and a central locking system, accessible from either side of the table, to secure table position.

Additional accessories are available, including foot controller, vascular scanning arm board, folding side rails, foot support, IV pole and paper dispenser.

SPECIFICATIONS:

Dimensions: 73" l x 30" w (185.4 x 76.2 cm); 35" w (88.9 cm)

with optional side rails

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Two primary sections: torso section features sonographer and patient cutouts with Fowler back positioning, and lower body section remains fixed

Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest

Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient

Mattress: Torso section 2" (5 cm) thick, lower body section 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) Individual locking swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 335 lb (152 kg)


Power: 115 VAC or 230 VAC

Warranty: One year

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2 No. 601.1-M90 and EN60601-1 standards.



More than 80% of sonographers are scanning in pain and 20% of these

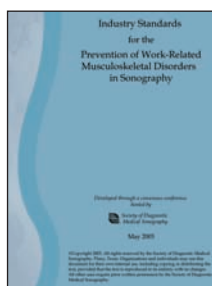

professionals eventually experience a career-ending injury. On average, within five years of entering the profession, sonographers experience pain while scanning.

For more information on SDMS Industry Standards for the Prevention of Work-Related Musculoskeletal Disorders in Sonography, visit our website at www.biodes.com

058-700 Table, Echocardiography, Echo Pro, 115

VAC

Includes Hand Controller.

*All options available separately**

058-705 Table, Echocardiography, Echo Pro, 230 VAC

Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-632 Side Rails, Folding

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-640 Foot Support

058-610 Dispenser, Paper

with Cutter Strap

Replacement:

056-612 Table Paper, Hygienic, 12 rolls

18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*

U.S. Patent No. 7,103,932 B1